




# RhythmReview guided demo

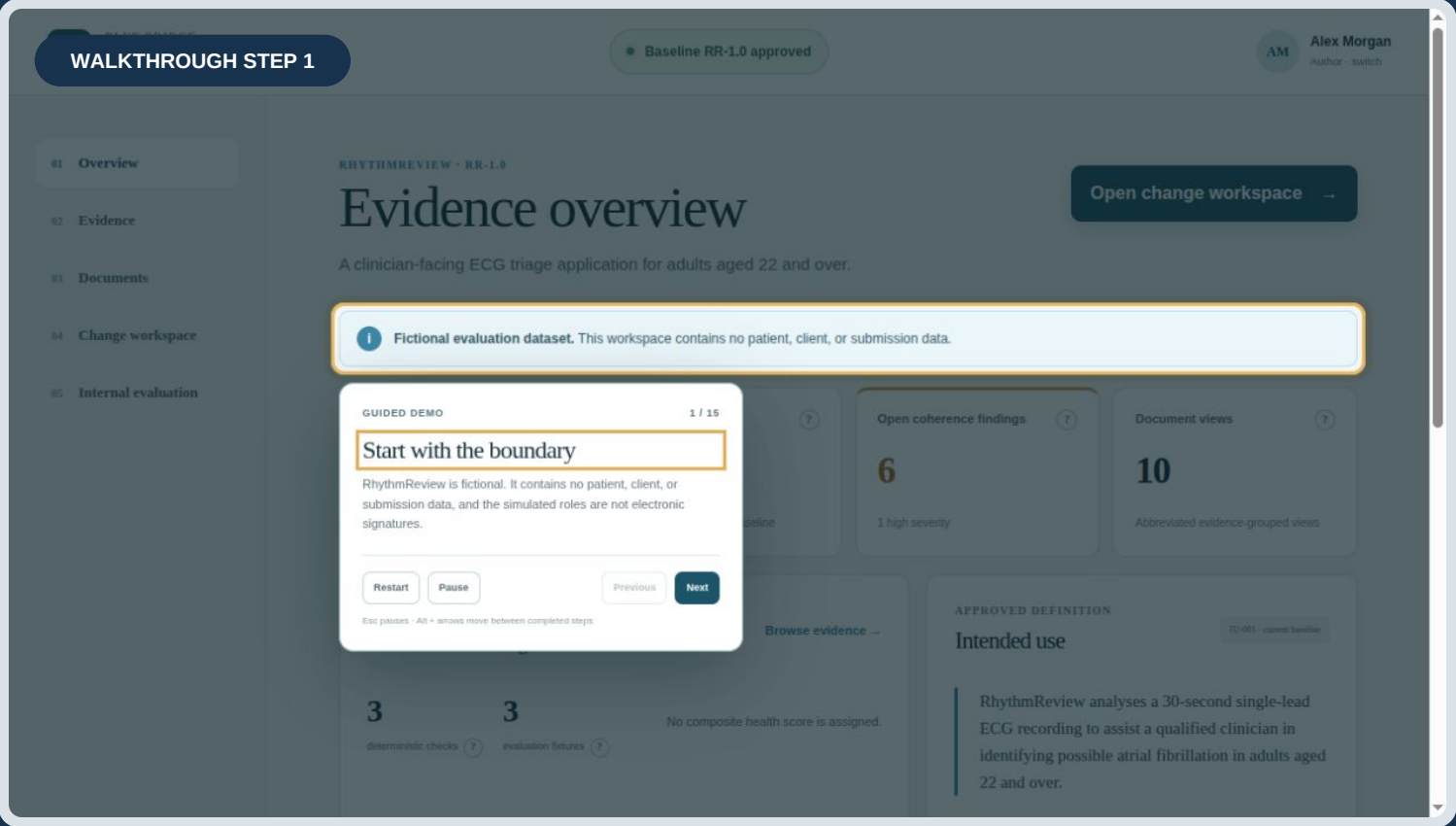

## Presentable for one purpose

A 10–12 minute internal review of workflow mechanics, human decision points, and immutable

**PROTOTYPE BOUNDARY**

Fictional data. Decision support only. Simulated roles are not electronic signatures.

Use the persistent “Start guided walkthrough” control.

The tour never performs a controlled action silently.

Reset the workspace before every presentation.

# Checks and evidence direction

2 minutes

The screenshot displays the RhythmReview interface during a guided demo. At the top, a 'WALKTHROUGH STEP 2' banner is visible. The main dashboard shows four key metrics: 72 (10 controlled types), 118 (in the active approved baseline), 6 (1 high severity), and 10 (Abbreviated evidence-grouped views). A central 'EXPLAINABLE CHECKS' panel highlights 'Coherence findings' with 3 deterministic checks and 3 evaluation fixtures, noting that no composite health score is assigned. A 'GUIDED DEMO' overlay titled 'Separate checks from fixtures' explains that three findings come from deterministic rules and three are seeded evaluation fixtures, avoiding a composite compliance score. The demo includes 'Restart', 'Pause', 'Previous', and 'Next' buttons. At the bottom, a 'RERUNNABLE CHECKS' section shows 'Approved baseline finding dispositions' and a 'Rerun checks' button.

## SAY THIS

Three findings are deterministic rules. Three are evaluation fixtures. Do not combine them into a compliance score.

- 2 Point to the two groups and explain why their provenance is visibly different.
- 3 Open REQ-004 next. Relationship labels report stored source and target direction only.

## QA/RA QUESTION

Which controlled dependency meanings should each relationship type carry?

# Direction and immutable snapshots

2 minutes

**WALKTHROUGH STEP 3**

Baseline RR-1.0 approved

**Clinical workflow study**  
Simulated-use study assessing result interpretation, acknowledgement, and referral decisions by clinicians.

**GUIDED DEMO** 3 / 15

**Read the stored relationship direction**

REQ-004 groups links by the saved source and target fields. Incoming and outgoing are factual storage directions; QA and RA still need to ratify dependency semantics.

Restart Pause Previous Next

Exit guides - Alt + arrow move between completed steps

**Engineering** High US + EU

**RATIONALE**  
Defines controlled requirement evidence for RhythmReview.

**INCOMING RELATIONSHIPS - THIS ITEM IS THE TARGET** 3

DES-004 — implements — REQ-004 Open DES-004

TEST-004 — verifies — REQ-004 Open TEST-004

DES-007 — may affect — REQ-004 Open DES-007

**OUTGOING RELATIONSHIPS - THIS ITEM IS THE SOURCE** 2

REQ-004 — refines — U7N-004 Open U7N-004

REQ-004 — may affect — HAZ-003 Open HAZ-003

**VERSION HISTORY**

v1.0 Approved Jamie Chen - QA reviewer

**Time to result**  
For supported recordings, the application returns a result within 30 seconds after upload completes.

**CLM-004 - Product claim**

**Clinician decision support**

**PROTOTYPE BOUNDARY**  
Decision support only. AI and replay fixtures cannot approve evidence or baselines.

Resume guided walkthrough 01-12-2024 - QA/RA

**WALKTHROUGH STEP 4**

Baseline RR-1.0 approved

**snapshot**

Each abbreviated document view records its baseline, source versions, approver, and snapshot ID. These views are not submission-ready documents.

Restart Pause Previous Next

Exit guides - Alt + arrow move between completed steps

**Controlled product purpose, users, population, claims, and exclusions.**

| Baseline status | Baseline | Approved by              | Snapshot               |
|-----------------|----------|--------------------------|------------------------|
| approved        | RR-1.0   | Jamie Chen - QA reviewer | SNAP-DOC-001-BL-RR-1.0 |

**01 Product claim**

CLM-001 - v1.0 RhythmReview identifies ECG recordings that may show atrial fibrillation for clinician review.

**Possible AF indication**

CLM-002 - v1.0 The application reports Possible AF. No AF detected, or Unreadable recording.

**Three controlled outcomes**

CLM-003 - v1.0 For supported recordings, the application returns a result within 30 seconds after upload completes.

**Time to result**

CLM-004 - v1.0 The result informs but does not replace clinical judgement or diagnostic confirmation.

**Clinician decision support**

**02 Intended use**

IU-001 - v1.0 RhythmReview analyses a 30-second single-lead ECG recording to assist a qualified clinician in identifying possible atrial fibrillation in adults aged 22 and over. It is not intended for continuous monitoring, emergency use, or standalone diagnosis.

## INCOMING / OUTGOING


Incoming means the current item is the stored target. Outgoing means it is the stored source. Neither label claims upstream or downstream dependency semantics.

## SNAPSHOT PROOF

Show baseline ID, snapshot ID, source versions, and approver. These are abbreviated controlled views, not submission-ready documents.

# Create the timing replay

2 minutes

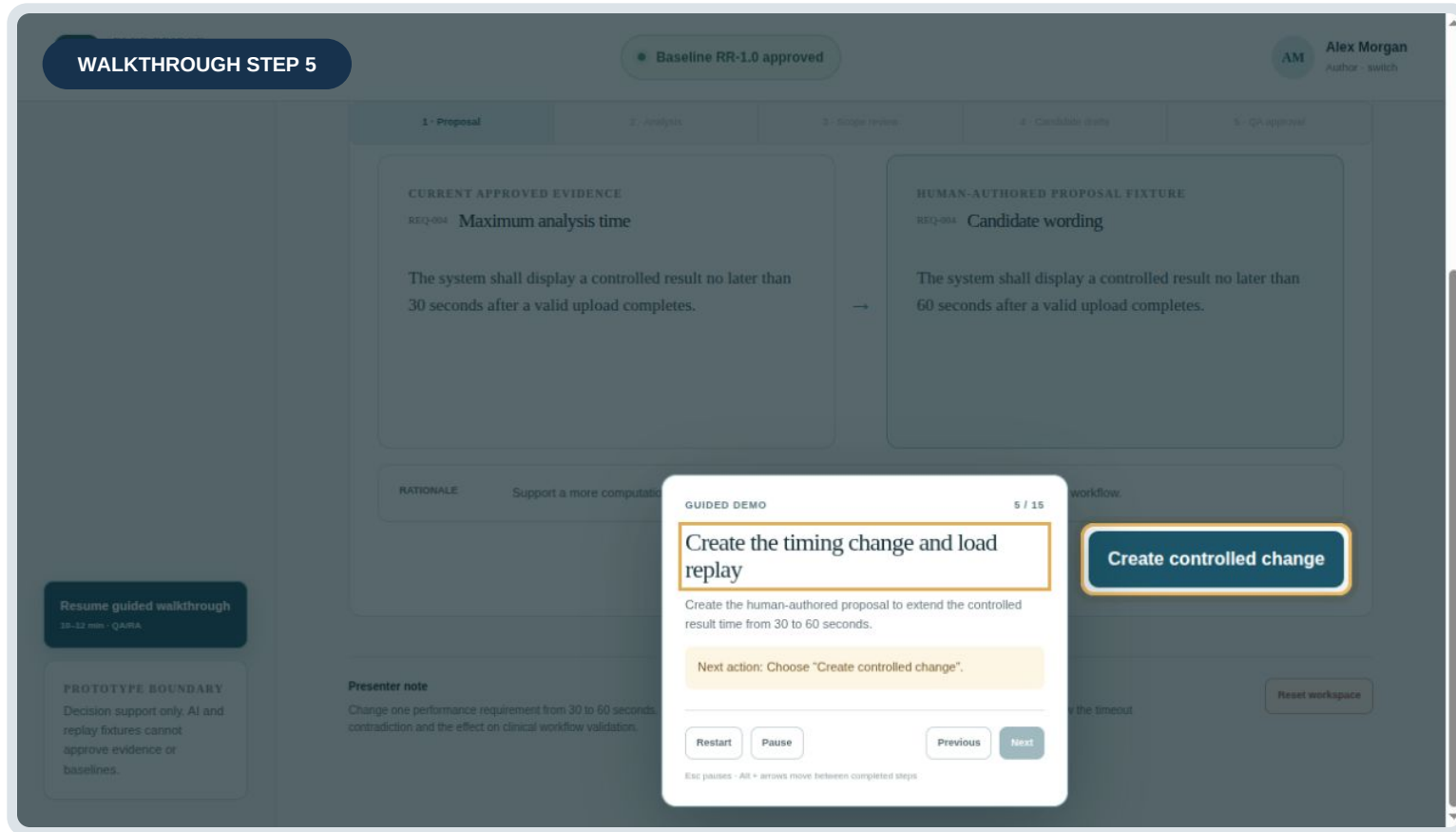

- 5 Create the human-authored proposal: extend the controlled result time from 30 to 60 seconds.
- 6 Load replay data. Compare deterministic graph paths with unlinked semantic candidates.
- 7 Switch to Jamie Chen. The role change remains an explicit presenter action.

## REPLAY BOUNDARY

Fixtures demonstrate workflow mechanics. They are not live model output and cannot approve, waive, or create controlled content.

# Make and retain human decisions

3 minutes

The screenshot displays the RhythmReview Phase 1.5 guided demo interface. A dark blue header bar at the top contains the text "WALKTHROUGH STEP 12" on the left and "GUIDED DEMO 12 / 15" on the right. The main content area is divided into two sections. The top section, titled "Author the candidate and rerun checks", includes a sub-header "Rerun the projected candidate checks after the author edit so the displayed findings are current." and a row of buttons: "Restart", "Pause", "Previous", and "Next". Below this is a section titled "RERUNNABLE CHECKS" with a sub-header "Projected candidate findings". It shows a summary "6 failing - 0 waived - Run 08/09/2026, 17:54:31" and a table of findings. The bottom section is titled "AUDIT HISTORY" and shows "10 recorded events".

**WALKTHROUGH STEP 12**

**GUIDED DEMO 12 / 15**

**Author the candidate and rerun checks**

Rerun the projected candidate checks after the author edit so the displayed findings are current.

Restart Pause Previous Next

Esc pauses - Alt + arrows move between completed steps

**RERUNNABLE CHECKS**

**Projected candidate findings**

6 failing - 0 waived - Run 08/09/2026, 17:54:31

|         |                                                       |      |
|---------|-------------------------------------------------------|------|
| CLM-006 | Label overstates the claim                            | OPEN |
| CHP-004 | Architecture component absent from inventory          | OPEN |
| REV-003 | Evidence review date missing                          | OPEN |
| TRC-014 | Risk control missing verification link                | OPEN |
| VAL-022 | Inconsistent timing value                             | OPEN |
| VER-009 | Superseded evidence appears in current document views | OPEN |


**AUDIT HISTORY**

10 recorded events


- 8 Accept REQ-004.
- 9 Reject TEST-007, then revise it to accepted. Open audit history to show both decisions.
- 10 Edit UN-004 from review to update.
- 11 Use the guided helper only after those examples. It fills pending decisions with explicit fixture provenance.
- 12 Switch to the author, edit a candidate, and rerun the projected candidate checks.

# Approve, audit, and compare

2 minutes

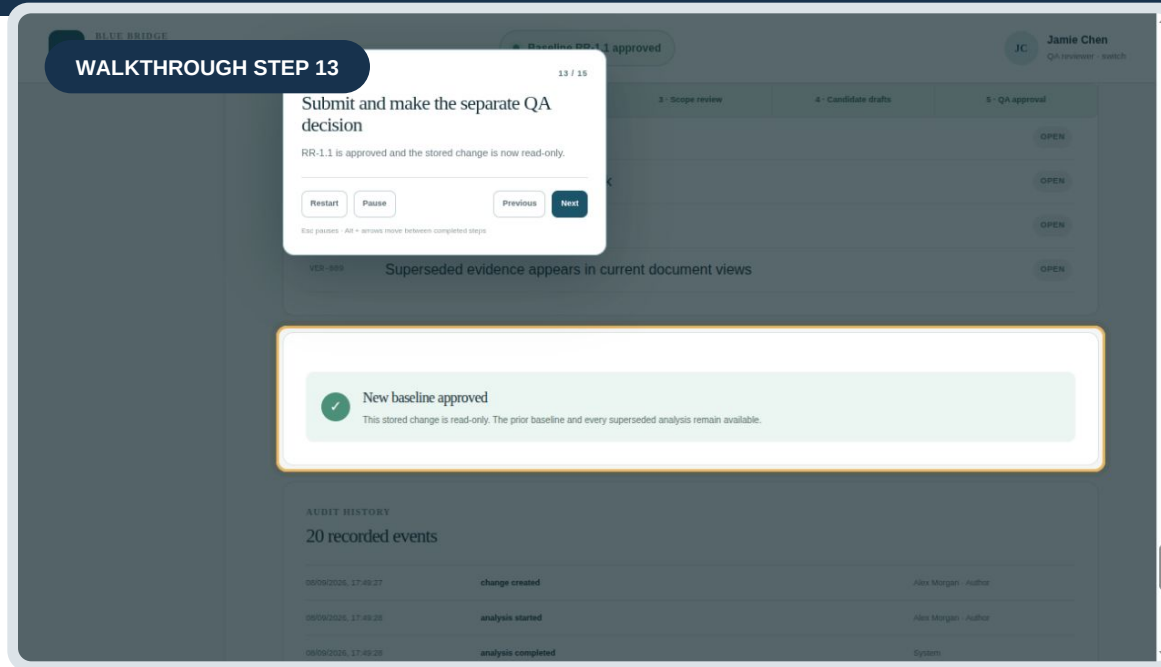

- 13 Submit as the author, switch to QA, then approve RR-1.1 as a separate human decision.

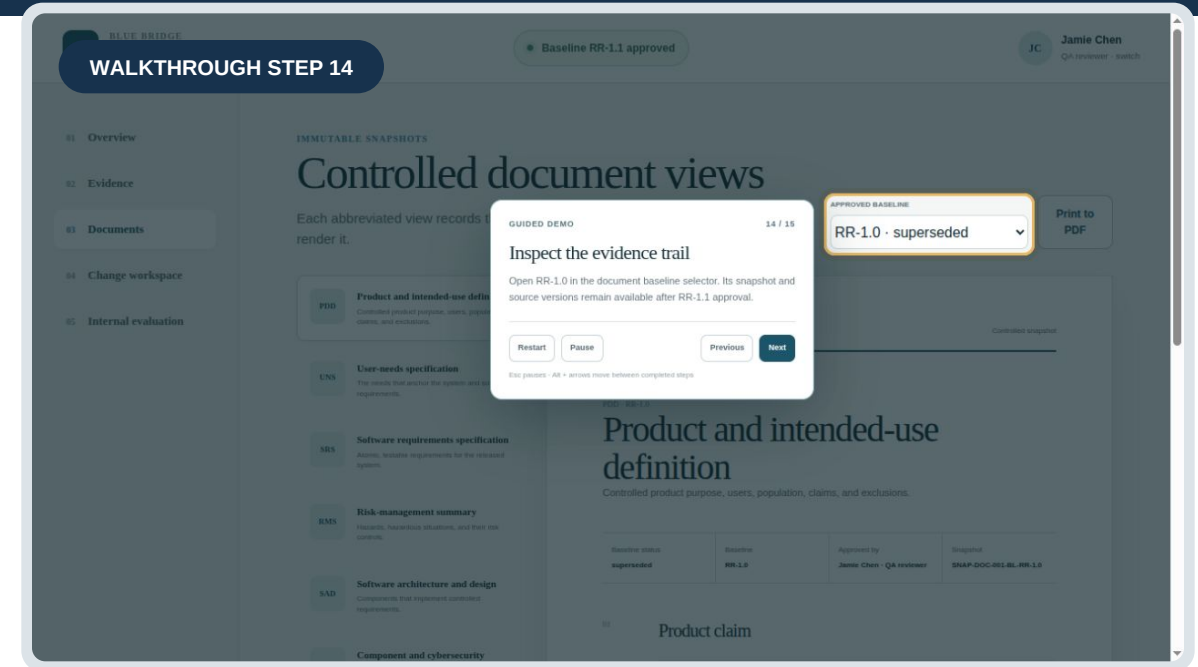

- 14 Expand an audit event, then select RR-1.0 to prove the earlier immutable snapshot is retained.

## RESET AFTER THE COMPARISON

Choose Reset workspace. The demonstration change is removed and RR-1.0 becomes current again.

# Evaluation is not performance evidence

Keep outside the main tour

**EVALUATION APPENDIX**

Baseline RR-1.0 approved

AM Alex Morgan  
Author · switch

01 Overview  
02 Evidence  
03 Documents  
04 Change workspace  
05 Internal evaluation

INTERNAL VALIDATION WORKFLOW

## Evaluation dashboard

This page is separate from the customer compliance workflow. It compares a stored run with a locked provisional human answer key.

SCENARIO: 01 - Clarify result wording

STORED STRUCTURED RUN: Latest completed run

**Not model-performance evidence.**  
Replay suggestions were authored from the same provisional answer key used here. Use this page to inspect calculation mechanics only until independent QA/RA review is complete.

**No completed structured run**  
Run a guided replay or live analysis, then record QA decisions. This page does not substitute seeded scores.

ILLUSTRATIVE COMPARISON FIXTURES

### Manual and document-chat values

Not calculated from this run

WORKFLOW RECALL PRECISION REVIEW TIME

## USE WITH CAUTION

Replay suggestions were authored from the same provisional answer key used for scoring. This can show calculation mechanics only.

## Open answer-key review

- ? CLM-003 is absent from scenario 2.
- ? LBL-004 is absent from scenario 3.
- ? REQ-001 has a questionable expected action.

**Independent QA/RA review is required before any performance claim.**

# Reset, review, then begin Phase 2

## Before handing over

- 1 Reset the workspace and confirm RR-1.0 is current.
- 2 Rehearse the 15-step path at presentation size.
- 3 State the prototype and replay boundaries before showing findings.
- 4 Keep evaluation in the appendix until the answer key is independently reviewed.
- 5 Record QA/RA decisions on relationship semantics and answer-key questions.

## Phase 2 roadmap

- 01 Traceability coverage matrix
- 02 Risk workflow
- 03 Verification workflow
- 04 Release readiness
- 05 Relationship editing

### AFTER PHASE 2

Add a fixture-backed Jira adapter. Live Jira remains outside Phase 1.5.